

Cleaning afterwards

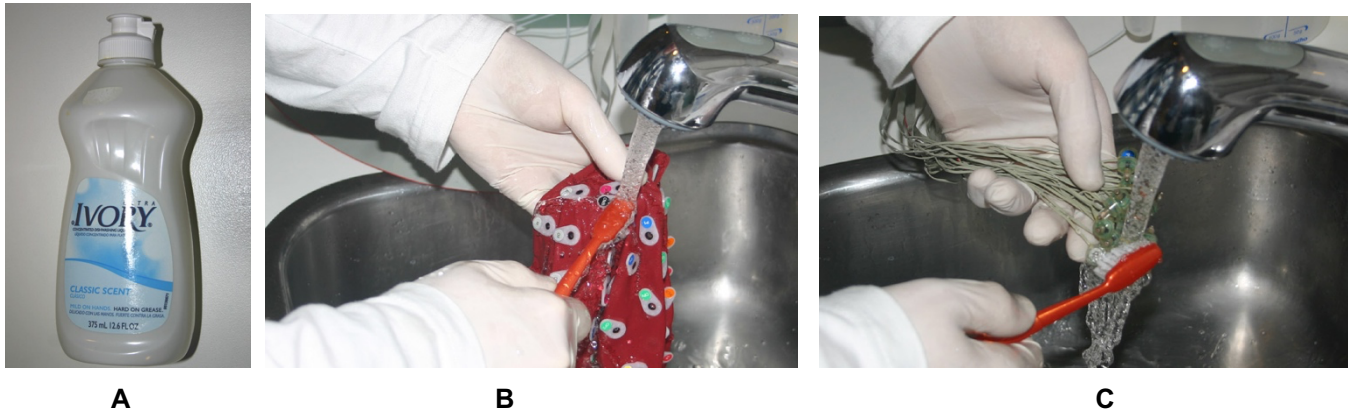


Figure 11: A) Ivory detergent is a regular dish washing soap that is recommended for initial cleaning of electrodes after use, because it is not abrasive (i.e., it does not scratch off the sensitive surface of the electrode). B) A wooden stick comes handy to remove the gel from all corners of the cap. Use the soap in combination with a soft toothbrush or dental pick to further remove residues of electrode gel from the cap and the electrodes before disinfection. C) Do not let the electrodes get in contact with metal. You can use a plastic box with some detergent to save water and to keep electrodes away from the metal in the kitchen sink.



Figure 12: A) In the kitchen, you can find dedicated containers for cleaning and disposing materials. Used needles should be put in a needle container. B) Sekusept Plus is the disinfectant that we use for cleaning electrodes and electrode caps after use. Make sure you have read the instructions before you use it. Prepare 500 ml of new solution of Sekusept in water in the round plastic container for one test day, or prepare a solution for multiple days in the black container and pour a day's portion into the round plastic container. Unused solution can be held for 14 days. C) If you use Sekusept in a 4% solution, keep the rinsed electrodes and electrode cap in the solution for 15 minutes – not longer, not shorter. The Velcro-straps should not be disinfected. After removing the electrode and cap from the solution, keep it under running water for some time to remove the aggressive chemicals that might otherwise come in contact with the next participant's hair or skin. Protect your own skin against the aggressive chemicals by wearing (new) gloves and protective glasses.